CUSTOMISED CLINICAL IMAGE IES - MULTILEVEL




CASTELLINI IES FILTERS (IMAGE ENHANCEMENT SYSTEM)


Latest-generation X-VS image processing software aims to improve diagnostics. With excellent image resolution and a user-friendly software interface, X-VS makes reading intraoral images simpler and easier. The new CASTELLINI iES (Image Enhancement System) filters are the result of research that focuses on dentists' real needs. By using proprietary algorithms optimised for the X-VS sensor, dentists can simultaneously capture, display and share a set of Multi Level images (up to 5). Each image is the result of a different improvement designed to highlight various anatomical details with different levels of sharpness. Dentists can customise image contrast to suit their diagnostic or visual preferences, allowing improved diagnosis. Afterwards, favourite settings can be rendered automatic.


Equipped with advanced iRYS software, the X-VS now offers the versatile MultiLEVEL function, which lets users pre-set the image processing filters displayed in Multi Level. Users can select which filters to use from among the pre-set families and define any further customisations, all from the iRYS image display window. This provides individual dentists with a customised comfort zone for every appointment.


41.9 x 30.4 External dimensions (mm) 38.9 x 24.9 Thickness (mm) 5.3 1500 x 1000 1700 x 1300 Pixel size (µm) 20 Maximum resolution (lp/mm) Grey levels depth 14-bit acquisition - 16384 maximum grey levels Csl (Cesium lodide) with micro-columnar structure Scintillator technology Direct exposure protection FOP (Fibre Optics Plate) Protection rating IP 67 (Guaranteed against liquid or dust infiltration) Compatibility with X-ray generators Any AC or DC technology X-ray generator with values in the 60 - 70 kV and 1-8 mA range and precision control of exposure times iCapture with dedicated filters for third party software Image capture software (for PC) Image management software (for PC) iRYS (as per ISDP®10003:2020 in compliance with EN ISO/IEC17065:2012 - certificate number 2019003109-2) and iPad iRYS viewer app (free) DICOM 3.0, TWAIN, VDDS Supported protocols IHE compliant (Print; Storage Commitment, SR document; WorkList; MPPS; Query/Retrieve) DICOM nodes Minimum system requisites Microsoft® Windows® 10 Professional 64 bit Supported operating systems 1280x1024; 1344 x768 or greater, 16 million colours Display settings USB 2.0 or subsequent Power supply 5 V DC. 500 mA (via USB) Constant potential, microprocessor-controlled Working frequency 145 - 230 KHz with self-adjustment (typically 175 KHz) Focal spot 0.4 mm (IEC 336) 2 mm @ 60 kV / 2 mm @ 65 kV / 2 mm @ 70 kV [Total filtration Anode current 4 / 8 mA 60 / 65 / 70 kV (*) Voltage at X-ray tube 0.020 - 1.000 seconds, R'10 and R'20 scale Exposure times Source-skin distance 20 and 30 cm Irradiated field Ø 60 mm and Ø 55 mm (with round cone) 35 x 45 mm (with rectangular cone for size 2 sensors) Additional collimators 31 x 41 mm and 22 x 35 mm, for size 1 and size 0 sensors 50/60 Hz, 115-120 V AC ±10% or 230-240 V AC ±10% Power supply Duty Cycle Continuous operation with self-adjustment up to 1s/90s tota Arms (for Standard version only) Available in 3 lengths: 40 cm - 60 cm - 90 cm Max. arm extension 230 cm, from wall Standard (wall mounted) or Mobile (on portable cart) **(€**

(*) values depend on the country where the product is marketed.

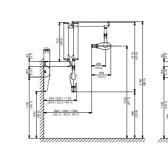

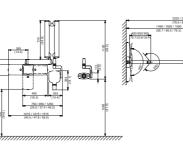

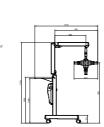

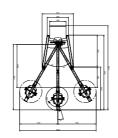



BU Medical Equipment Sede Legale Ed Amministrativa Headquarters

CEFLA s.c. - Via Selice Provinciale, 23/a - 40026 Imola (BO) - Italy Tel +39 0542 653111 Fax +39 0542 653344

Stabilimento

CEFLA s.c. - Via Bicocca, 14/C - 40026 Imola (BO) - Italy Tel. +39 0542 653441

Fax +39 0542 653601

CASTELLINI.COM

IMAGING RX DC x-vs

THE PERFECT DIMENSION OF IMAGING






RX DC X-VS IMAGING THE PERFECT DIMENSION OF IMAGING

IMMEDIATE DIAGNOSIS, **EXCELLENT RESULTS**



The RX DC is a versatile, user-friendly X-ray unit, capable of producing high quality imaging thanks to cutting-edge technology. It also maximises working comfort while ensuring low X-ray doses for the patient and maintaining ultra-high performance. The RX DC, in fact, uses a constant potential high frequency (DC) generator and a very small focal spot (0.4 mm) that provides sharp, detailed images. Automatic exposure parameter modulation and accurate power selection make the RX DC the perfect X-ray unit whatever the situation, personalised according to the patient's build and the specific region of investigation. The RX DC is the perfect X-ray unit, capable of combining high quality imaging and a versatile, ergonomic design with low patient doses.

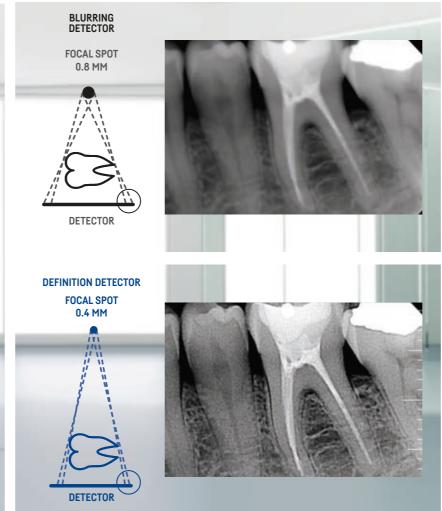

MINIMUM DOSE, HIGH DEFINITION, HIGH FREQUENCY RX DC X-RAY UNIT

The RX DC increases X-ray parallelism: thanks to the incorporated collimator the RX DC can achieve a source-to-skin gap of 30 cm. Sharp edges, clear images and precise detail, all while keeping patient doses low. Ergonomic design, offers simplicity and reliability thanks to solid arms with an integrated self-balancing system. Arm and tube head positioning is comfortable and stable. The protractor with graduated scale allows easy positioning.

MAXIMUM VERSATILITY AND MOBILITY

The RX DC is now even more practical and versatile, wall-mounted with 6 variable positions - 3 extensions available in the following lengths: 40, 60 and 90 cm - or, thanks to the cart, it can easily be moved around for use in different areas of the surgery. One X-ray unit for all your needs.





COLLIMATOR CONES

With the built-in collimator cone, a focusskin distance of 30 cm can be reached. The optional rectangular cone further reduces the body area exposed to X-rays.



USER-FRIENDLY CONTROL

Simple, user-friendly handheld unit, can be used to select the most suitable programme and ensure perfect X-ray acquisition. Moreover, thanks to the fast Dynamic Duty-Cycle, it is possible to keep tube temperature under control and check - in real time - the exact administered dose.


Innovative ergonomics, direct USB plug-and-play connection, high definition and immediate results make the X-VS with X-ViSUS HD technology the most advanced and suitable sensor for your surgery. Simplicity of use and image acquisition - combined with advanced real-time digital technology - improve quality of work. Impact-resistant, dust-resistant and certified IP67 (water-resistant), the X-VS - X-ViSUS HD technology can be used with all X-ray systems. The X-VS - X-ViSUS HD technology uses iRYS, the all-in-one software for diagnostics, communication and management of intraoral imaging which allows flawless storage, processing and printing of images in perfect synergy with any other devices already in the surgery.





MADE-TO-MEASURE SENSOR

of the patient's oral cavity. Outstanding positioning comfort thanks to rounded and a maximised active area.

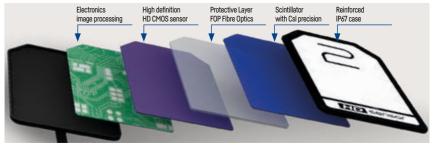

LATEST-GENERATION HD SENSOR

X-VS - X-ViSUS HD technology allows users 4-layer sensor with an additional protective layer to provide sharp, high-contrast images. to choose between two sensor sizes; these Caesium lodide (CsI) scintillator made up of column-like micro-structures that preserve make it more adaptable to the dimensions image quality; it first intercepts the X-ray beam and converts it into visible light. The Fibre Optics Plate (a layer of fibre optics) collimates the radiation onto the sensor and protects it against direct X-ray penetration. The high definition (CMOS HD) acquisition device and oncorners; high performance, compact design board electronics convert the light into a digital image with 16,384 grey levels.

RX DC X-VS IMAGING THE PERFECT DIMENSION OF IMAGING RX DC X-VS IMAGING THE PERFECT DIMENSION OF IMAGING