
MILLSEAL ROTECTION FOR YOUR THERMAL SEALERS INSTRUMENTS

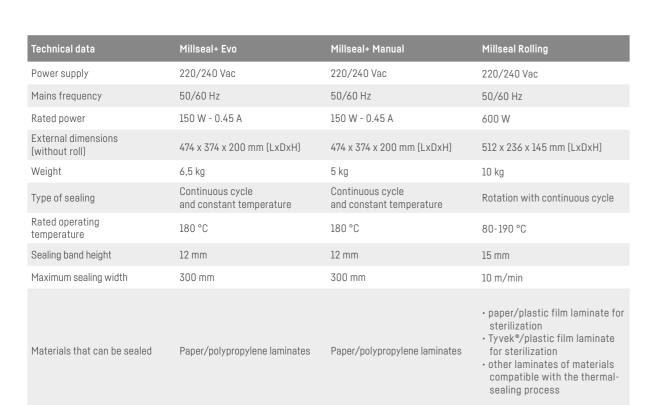

200 mm

Castellini sealers conform to the applicable standards:

Maximum roll diameter

200 mm



MILLSEAL THERMAL SEALERS PROTECTION FOR YOUR INSTRUMENTS





Application of CE 2006/42, CE 2006/95, 2004/108 Directives.
Reference standards EN61010-1:2001, EN61326-1:2006

A MORE STREAMLINED WORKFLOW

MILLSEAL THERMAL SEALERS

Castellini thermal sealers offer an efficient and reliable solution for dental surgeries. The Millseal + Manual and Millseal + Evo models, with built-in cutters, are practical and ergonomic. The Millseal + Evo version, in particular, has user-friendly automatic controls to manage the cutting and automatic first sealing of large numbers of pouches. The Millseal Rolling thermal sealer completes the range: it features an automatic feed system to optimize the speed of sealing operations and save considerable time for the sterilization room staff.



MILLSEAL+ EVO AUTOMATIC THERMAL SEALER

Pouch length and number indicated on the display. Two-way roller slide system. Cut and first pouch seal key. Definitive pouch seal key. Program keys allow 2 different pouch lengths to be saved and sealing/cutting tasks to be executed automatically.

A PROGRAMME FOR EVERY NEED

Millseal+ Evo lets users set two programmes with different pouch lengths and quantities. Once a programme is selected, the thermal sealer automatically prepares the pouches. Alternatively, the operator can package the instruments manually

MILLSEAL THERMAL SEALERS

in a more traditional manner.





KEYPAD BAG HOLDER

PROTECTION FOR YOUR INSTRUMENTS



MILLSEAL+ MANUAL THERMAL SEALER

The selection and cutting of the pouch length are performed manually. To seal the pouch containing the material ready for sterilization, simply press the welding lever until the acoustic signal indicates completion of the sealing. A positioning systems holds the paper in place to prevent rewinding of the roll after cutting and sealing.

MILLSEAL ROLLING THERMAL SEALER

Roller-type thermal sealer with automatic feed system. This model offers medical-standard quality and is ideal for speeding up and optimising the sealing of pouches that are pre-cut and sealed on 3 sides. The 15 mm sealing band ensures the long-lasting sterility of the packaged instruments. Temperature and sealing pressure adjustment provides effective control of the sealing process, for consistent high performance. The large display allows operating parameters to be set and viewed in an intuitive, user-friendly manner. The practicality of the product is enhanced by the various accessories (roller, work surface, cutter).









STAINLESS STEEL WORK SURFACE LARGE DISPLAY

ROLLER